

Smooth Muscle Cells

Instruction Manual

Product	Size	Catalog Number
Human Aortic Smooth Muscle Cells (HAoSMC)	500,000 cryopreserved cells	C-12533
	500,000 proliferating cells	C-12532
Human Coronary Artery Smooth Muscle Cells (HCASMC)	500,000 cryopreserved cells	C-12511
	500,000 proliferating cells	C-12512
Human Pulmonary Artery Smooth Muscle Cells (HPASMC)	500,000 cryopreserved cells	C-12521
	500,000 proliferating cells	C-12522
Human Umbilical Artery Smooth Muscle Cells (HUASMC)	500,000 cryopreserved cells	C-12500
	500,000 proliferating cells	C-12550
Human Tracheal Smooth Muscle Cells (HTSMC)	500,000 cryopreserved cells	C-12565
	500,000 proliferating cells	C-12566
Human Bronchial Smooth Muscle Cells (HBSMC)	500,000 cryopreserved cells	C-12561
	500,000 proliferating cells	C-12562
Human Bladder Smooth Muscle Cells (HBdSMC)	500,000 cryopreserved cells	C-12571
	500,000 proliferating cells	C-12572
Human Prostate Smooth Muscle Cells (HPrSMC)	500,000 cryopreserved cells	C-12573
	500,000 proliferating cells	C-12574
Human Uterine Smooth Muscle Cells (HUTSMC)	500,000 cryopreserved cells	C-12575
	500,000 proliferating cells	C-12576

Product Description

Smooth muscle tissue is found in the *tunica media* layer of large and small blood vessels and in the walls of hollow organs like the bladder and the uterus. Smooth Muscle Cells (SMC) possess all the same types of filaments, but, depending on the tissue of origin, differ significantly in mechanical and physiological properties. Therefore, PromoCell offers a range of Human Smooth Muscle Cells produced at PromoCell's cell culture facility from normal human tissue of different origins. The cells are isolated from the aorta, the coronary artery, the pulmonary artery, the umbilical artery, the trachea, the bronchi, the bladder, and the uterus.

Shortly after isolation, all PromoCell Human Smooth Muscle Cells are cryopreserved at passage 2 (P2) using PromoCell's proprietary, serum-free freezing medium, Cryo-SFM. Each cryo vial contains more than 500,000 viable cells after thawing.

Proliferating cell cultures are made from 500,000 cryopreserved cells that have been thawed and cultured for three days at PromoCell.

Quality Control

Rigid quality control tests are performed for each lot of PromoCell Human Smooth Muscle Cells.

They are tested for cell morphology, adherence rate, and cell viability. Immunohistochemical tests for cell-type specific markers, e.g. smooth muscle α -Actin (see page 5), are carried out for each lot. Growth performance is tested through multiple passages up to 15 population doublings (PD) under culture conditions without antibiotics and antimycotics.

In addition, all cells have been tested for the absence of HIV-1, HBV, HCV, and microbial contaminants (fungi, bacteria, and mycoplasma).

A detailed certificate of analysis (CoA) for each lot can be downloaded at: www.promocell.com/coa

Intended Use

PromoCell Human Smooth Muscle Cells are for *in vitro* research use only and not for diagnostic or therapeutic procedures.

Warning

Although tested negative for HIV-1, HBV, and HCV, the cells - like all products of human origin - should be handled as potentially infectious. No test procedure can completely guarantee the absence of infectious agents.

Follow appropriate safety precautions!

After delivery, start immediately with the protocol for cryopreserved cells (see page 2) or the protocol for proliferating cells (see page 3).

*Start immediately after delivery.
Use aseptic techniques and a laminar flow bench.*

Protocol for Cryopreserved Cells

Straight after arrival, store the cryopreserved cells in liquid nitrogen, or seed them immediately.

Note: Storage at -80°C is not sufficient for cell preservation and causes irreversible cell damage.

1. Prepare the medium

Calculate the needed culture surface area according to the plating density (see page 5). Fill the appropriate volume of PromoCell Growth Medium (at least 9 ml per vial of cells) in cell culture vessels. Place the vessels in an incubator (37°C , 5% CO_2) for 30 minutes.




2. Thaw the cells

Remove the cryovial from the liquid nitrogen container and immediately place it on dry ice - even for short transportation. Submerge the vial into a water bath (37°C) and continuously agitate for 90 sec.



3. Disinfect the vial and seed the cells

Thoroughly rinse the cryovial with 70% ethanol to avoid microbial contamination. Then, wipe the vial with a tissue. Open the vial under a laminar flow bench and resuspend the cells by carefully pipetting up and down. Transfer the cells to a cell culture vessel with the prewarmed medium from step 1.

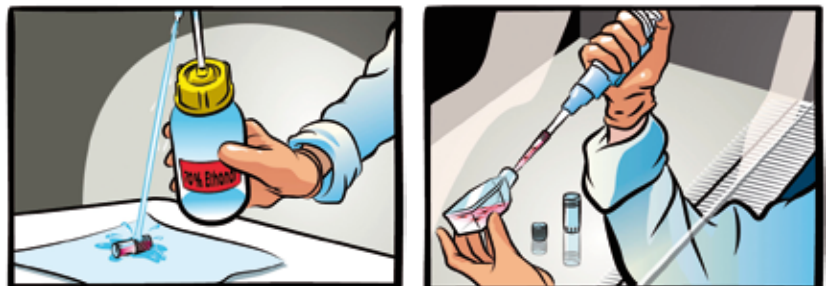

4. Incubate the cells

Place the vessel in an incubator (37°C , 5% CO_2) for cell attachment. Replace the medium after 16 – 24 hours. The cells should be subcultured, according to the subcultivation protocol (see page 4), once they have reached 70 – 90% confluency.



*Start immediately after delivery.
Use aseptic techniques and a laminar flow bench.*

Protocol for Proliferating Cells


1. Incubate the cells

Unpack the culture vessel, do not open the lid, and immediately place it in an incubator (37°C, 5% CO₂) for 3 hours.




2. Replace the transport medium

Carefully open the vessel, rinse the inner side of the lid with 70% ethanol, and let air dry. Aspirate the transport medium from the vessel. Add 10 ml of the appropriate PromoCell Cell Growth Medium.

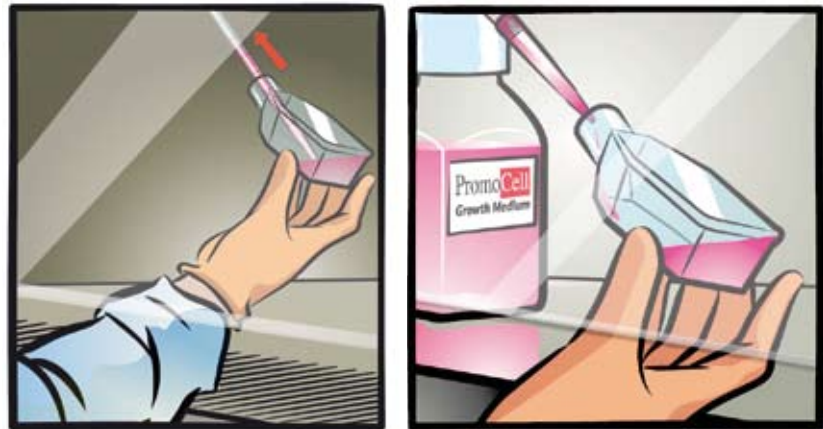

3. Check and incubate the cells

Check the cell density. Open the lid a half turn and place the vessel in an incubator (37°C, 5% CO₂). The cells should be subcultured, according to the subcultivation protocol (see page 4), once they have reached 70 – 90% confluency.

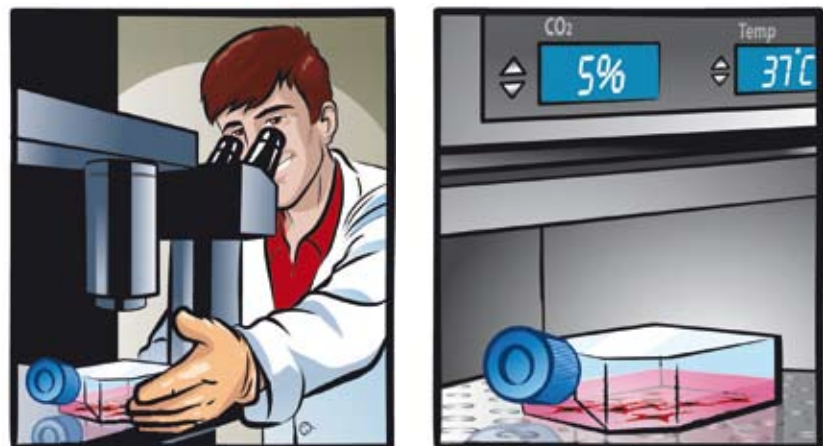

Use aseptic techniques and a laminar flow bench.

Subcultivation Protocol

1. Prepare the reagents and wash the cells

Place the PromoCell DetachKit at room temperature for at least 30 minutes to adjust the temperature of the reagents. Carefully aspirate the medium from the culture vessel. Add 100 μ l Hepes BSS Solution per cm^2 of vessel surface to wash the cells and agitate the vessel carefully for 15 seconds.



2. Detach the cells

Carefully aspirate the Hepes BSS from the culture vessel. Add 100 μ l Trypsin/EDTA Solution per cm^2 of vessel surface. Note: We recommend detaching the cells at room temperature. Close the vessel and examine the cells under a microscope. When the cells start to detach, gently tap the side of the vessel to loosen the remaining cells.




3. Neutralize the trypsin and harvest the cells

Add 100 μ l Trypsin Neutralization Solution per cm^2 of vessel surface and gently agitate. Carefully aspirate the cell suspension and transfer it to a centrifugation tube. Spin down the cells for 3 minutes at 220 x g.



4. Incubate the cells

Discard the supernatant (step 1), add 1 ml of the appropriate PromoCell Cell Growth Medium (step 2), and resuspend the cells by carefully pipetting up and down. Plate the cells according to the recommended seeding density in new cell culture vessels containing PromoCell Cell Growth Medium prewarmed to 37°C. Place the vessels in an incubator (37°C, 5% CO_2).

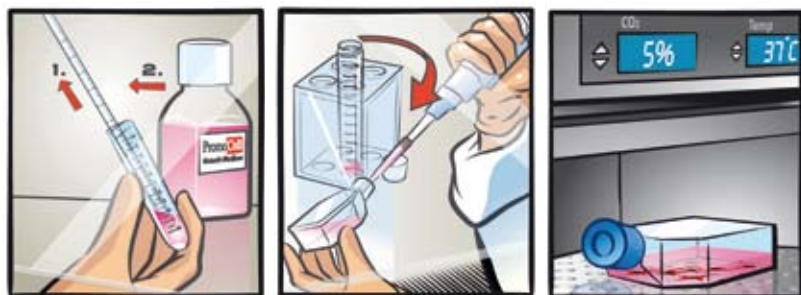

Specifications

Product	Recommended Culture Media*	Plating Density	Passage after Thawing	Marker	Population Doublings
Human Aortic Smooth Muscle Cells (HAoSMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15
Human Coronary Artery Smooth Muscle Cells (HCASMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15
Human Pulmonary Artery Smooth Muscle Cells (HPASMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15
Human Umbilical Artery Smooth Muscle Cells (HUASMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15
Human Tracheal Smooth Muscle Cells (HTSMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15
Human Bronchial Smooth Muscle Cells (HBSMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15
Human Bladder Smooth Muscle Cells (HBdSMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15
Human Prostate Smooth Muscle Cells (HPrSMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15
Human Uterine Smooth Muscle Cells (HUTSMC)	C-22062	7,500 - 10,000 cells per cm ²	P2	Smooth muscle α -Actin ⁺	> 15

*The catalog numbers in this table are for media in ready-to-use packaging.

Related Products

Product	Size	Catalog Number
Smooth Muscle Cell Growth Medium 2 (Ready-to-use)	500 ml	C-22062
Smooth Muscle Cell Growth Medium 2 Kit	500 ml	C-22162
Smooth Muscle Cell Basal Medium 2	500 ml	C-22262
Smooth Muscle Cell Basal Medium 2, phenol red-free	500 ml	C-22267
Smooth Muscle Cell Growth Medium 2 SupplementMix	for 500 ml	C-39267
Smooth Muscle Cell Growth Medium 2 SupplementPack	for 500 ml	C-39262
DetachKit	30 ml 125 ml 250 ml	C-41200 C-41210 C-41220
Cryo-SFM	30 ml 125 ml	C-29910 C-29912
HAoSMC Pellet	> 1 million cells per pellet	C-14053
HCASMC Pellet	> 1 million cells per pellet	C-14052
HPASMC Pellet	> 1 million cells per pellet	C-14051
HUASMC Pellet	> 1 million cells per pellet	C-14050
HTSMC Pellet	> 1 million cells per pellet	C-14056
HBSMC Pellet	> 1 million cells per pellet	C-14055
HBdSMC Pellet	> 1 million cells per pellet	C-14057
HPrSMC Pellet	> 1 million cells per pellet	C-14058
HUtSMC Pellet	> 1 million cells per pellet	C-14059
Human Pericytes from Placenta (hPC-PL)	500,000 cryopreserved cells 500,000 proliferating cells	C-12980 C-12981

PromoCell GmbH

Sickingenstr. 63/65
69126 Heidelberg
Germany

Email: info@promocell.com
www.promocell.com

North America

Phone: 1 – 866 – 251 – 2860 (toll free)
Fax: 1 – 866 – 827 – 9219 (toll free)

Deutschland

Telefon: 0800 – 776 66 23 (gebührenfrei)
Fax: 0800 – 100 83 06 (gebührenfrei)

France

Téléphone: 0800 90 93 32 (ligne verte)
Téléfax: 0800 90 27 36 (ligne verte)

United Kingdom

Phone: 0800 – 96 03 33 (toll free)
Fax: 0800 – 169 85 54 (toll free)

Other Countries

Phone: +49 6221 – 649 34 0
Fax: +49 6221 – 649 34 40